

# CodeLink™ 16 Assay Accessory Kit (310032)


---

Part #	Description	Qty Ea
600037	Centrifuge Drying Trays, 3 Chamber	2
600038	Counterbalances for Centrifuge Drying Trays	2
600039	16 Assay Sealing Tool	3

**Centrifuge Drying Trays (600037):** hold three slides with 16- Assay Chambers to allow for consistent processing and for drying of the arrays in the centrifuge.



**Counterbalances for Centrifuge Drying Tray (600038):** Required for drying by centrifugation if number of slides processed is not even, otherwise centrifuge will not balance. This item is in the same shape as a slide with a 16-assay chamber, and fits into one of the slots for a chambered slide to balance the Centrifuge Drying Trays on each side of the rotor.

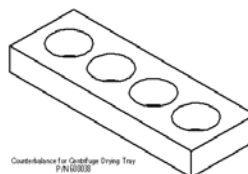

**16-Assay Sealing Tool (600039):** This tool is used to form a tight seal over the 16-Assay chamber for the incubation steps of the assay, particularly the hybridization, to ensure evaporation does not occur and alter the hybridization solution. This tool looks like a miniature rake.

